




Simplicity when you want it. Security when you need it.

Designed to provide a safe,
efficient, and reproducible
solution for anterior
cervical plating

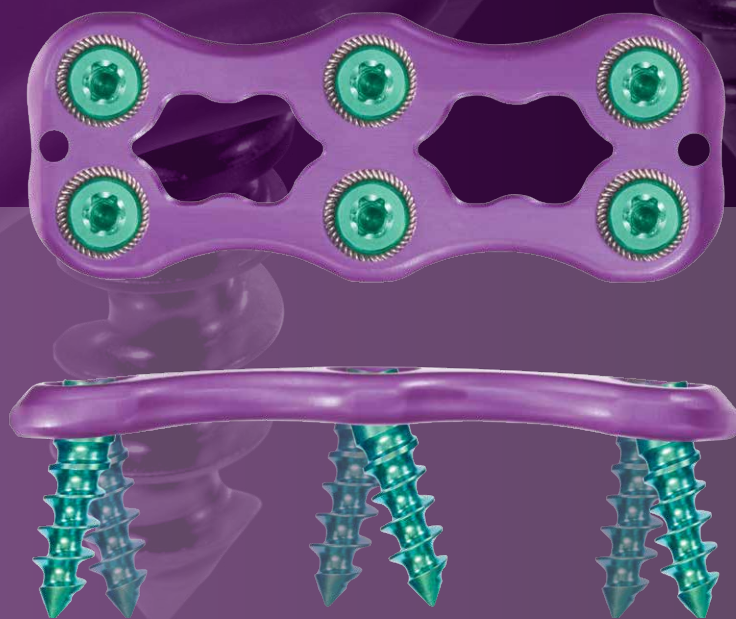

NuVasive® Helix ACP™ Overview

NuVasive Helix ACP

- NuVasive Helix ACP plate has a 2.5mm profile to reduce the chance of postoperative discomfort.
- 5-Level plates available via special order.
- Fixed- and Variable-Angle Screws for the perfect construct type, depending on patient needs.
- Variable-Angle Screws feature a 20° cone-of-angulation to accommodate surgeon needs.



NuVasive Helix ACP Canted Coil Lock (CCL)

- Consistently and securely blocks screws into the plate while still allowing the plate to be securely lagged to bone for a tight fit.
- Screws can be removed and repositioned easily while maintaining CCL integrity.
- Aggressive screws in a full complement of lengths and diameters accommodate a variety of surgical scenarios.



UNLOCKED

Ledge of screw is visible, and head of screw is proud to plate surface.




LOCKED

Ledge of screw is completely covered, and head of screw is seated below NuVasive Helix ACP Canted Coil Lock.

Balanced and Intuitive Instrumentation


- DTS Guides allow pre-drilling, tapping, and screw placement without changing insertion angle.

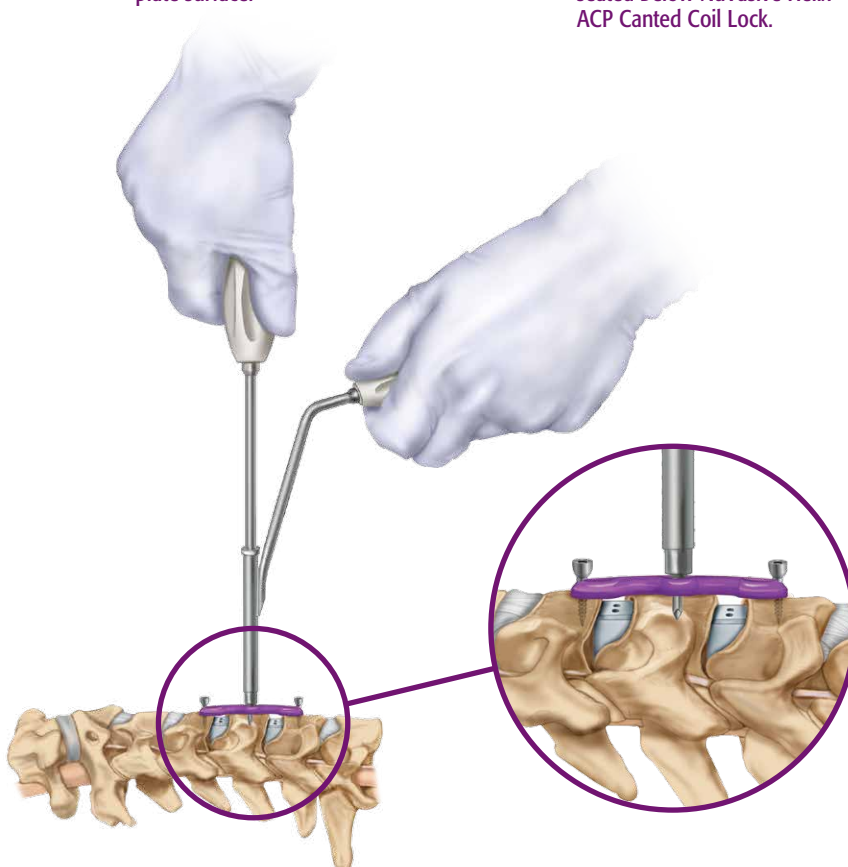

NuVasive® Helix ACP™ Overview

Low-Profile Leading Edge

- Plate sits flush against anterior cervical spine.
- Decreased “step” on leading edge of plate to reduce esophageal discomfort postoperatively.
- 1mm leading edge.

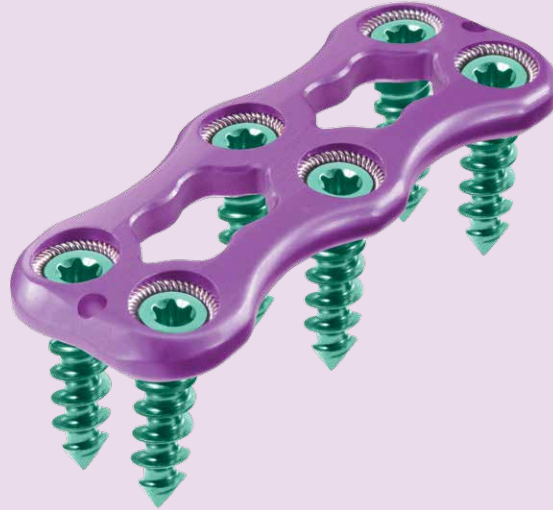

Narrow Plate Reduces Retraction During Placement

- Reduced retraction for easy placement.
- Less than 16mm at widest point.



Intraoperative Visualization

- Large graft windows allow for graft placement verification.



NUVASIVE® HELIX ACP™ INSTRUMENTS

DESCRIPTION	CATALOG #
13mm Drill (Sterile)	7730016
13mm Tap Bit	7730026
Screwdriver	7730041
Awl	7730062
Self-Centering Awl (Variable)	7730063
Self-Centering Awl (Fixed)	7730064
Screw Extractor	7730072
Drill Guide, Fixed	7730081
Drill Guide, Variable	7730091
Plate Bender	7730701
NuVasive Helix ACP Sterilization Tray	7730002
DTS Guide - Fixed, 0°	7730086
DTS Guide - Fixed, 7°	7730088
Plate Holder	7730043
Universal Handle	7730600
TEMPORARY FIXATION	
Temporary Fixation Pin (Temporary Tack)	7737110

NUVASIVE HELIX ACP IMPLANTS

DESCRIPTION	CATALOG #
SELF-DRILLING SCREWS (VARIABLE, 4.0)	
4.0 x 11mm Self-Drilling Variable Screw	7739111
4.0 x 13mm Self-Drilling Variable Screw	7739113
4.0 x 15mm Self-Drilling Variable Screw	7739115
SELF-TAPPING SCREWS (VARIABLE, 4.0)	
4.0 x 11mm Self-Tapping Variable Screw	7733111
4.0 x 13mm Self-Tapping Variable Screw	7739213
4.0 x 15mm Self-Tapping Variable Screw	7739215
4.0 x 17mm Self-Tapping Variable Screw	7739217
4.0 x 19mm Self-Tapping Variable Screw	7733119
SELF-TAPPING SCREWS (VARIABLE, 4.5)	
4.5 x 13mm Self-Tapping Variable Screw	7735113
4.5 x 15mm Self-Tapping Variable Screw	7735115
4.5 x 17mm Self-Tapping Variable Screw	7735117
4.5 x 19mm Self-Tapping Variable Screw	7735119
SELF-DRILLING SCREWS (FIXED, 4.0)	
4.0 x 11mm Self-Drilling Fixed Screw	7732111
4.0 x 13mm Self-Drilling Fixed Screw	7732113
4.0 x 15mm Self-Drilling Fixed Screw	7732115

NUVASIVE HELIX ACP IMPLANTS

DESCRIPTION	CATALOG #
SELF-TAPPING SCREWS (FIXED, 4.0)	
4.0 x 11mm Self-Tapping Fixed Screw	7734111
4.0 x 13mm Self-Tapping Fixed Screw	7734113
4.0 x 15mm Self-Tapping Fixed Screw	7734115
4.0 x 17mm Self-Tapping Fixed Screw	7734117
4.0 x 19mm Self-Tapping Fixed Screw	7734119
SELF-TAPPING SCREWS (FIXED, 4.5)	
4.5 x 13mm Self-Tapping Fixed Screw	7736113
4.5 x 15mm Self-Tapping Fixed Screw	7736115
4.5 x 17mm Self-Tapping Fixed Screw	7736117
PLATES	
1-Level NuVasive Helix ACP, 20mm	7738120
1-Level NuVasive Helix ACP, 22mm	7738122
1-Level NuVasive Helix ACP, 24mm	7738124
1-Level NuVasive Helix ACP, 26mm	7738126
1-Level NuVasive Helix ACP, 28mm	7738128
1-Level NuVasive Helix ACP, 30mm	7738130
1-Level NuVasive Helix ACP, 32mm	7738132
1-Level NuVasive Helix ACP, 34mm	7738134
2-Level NuVasive Helix ACP, 36mm	7738236
2-Level NuVasive Helix ACP, 38mm	7738238
2-Level NuVasive Helix ACP, 40mm	7738240
2-Level NuVasive Helix ACP, 42mm	7738242
2-Level NuVasive Helix ACP, 44mm	7738244
2-Level NuVasive Helix ACP, 46mm	7738246
2-Level NuVasive Helix ACP, 48mm	7738248
2-Level NuVasive Helix ACP, 50mm	7738250
2-Level NuVasive Helix ACP, 52mm	7738252
2-Level NuVasive Helix ACP, 54mm	7738254
3-Level NuVasive Helix ACP, 56mm	7738356
3-Level NuVasive Helix ACP, 58mm	7738358
3-Level NuVasive Helix ACP, 60mm	7738360
3-Level NuVasive Helix ACP, 62mm	7738362
3-Level NuVasive Helix ACP, 64mm	7738364
3-Level NuVasive Helix ACP, 66mm	7738366
3-Level NuVasive Helix ACP, 68mm	7738368
4-Level NuVasive Helix ACP, 70mm	7738470
4-Level NuVasive Helix ACP, 74mm	7738474
4-Level NuVasive Helix ACP, 78mm	7738478
4-Level NuVasive Helix ACP, 82mm	7738482
4-Level NuVasive Helix ACP, 86mm	7738486
4-Level NuVasive Helix ACP, 90mm	7738490

5-Level NuVasive Helix ACP plates are not part of the standard configuration, but may be ordered separately by contacting NuVasive Customer Service.



To order, please contact your NuVasive® Sales Consultant or Customer Service Representative today at:
NuVasive, Inc. 7475 Lusk Blvd., San Diego, CA 92121 USA • phone: 800-475-9131 fax: 800-475-9134
NuVasive UK Ltd. Suite B, Ground Floor, Caspian House, The Waterfront, Elstree, Herts WD6 3BS UK



phone: +44 (0) 208-238-7850 fax: +44 (0) 207-998-7818

www.nuvasive.com

